



# Quantification of meclofenamic acid in equine plasma by ultra-high performance liquid chromatography electrospray ionization tandem mass spectrometry (UHPLC-ESI/MS/MS)

Madison Chao, Shelby Carlson and Amari Spates

Faculty Mentor: Liguo Song

Department of Chemistry, Western Illinois University, Macomb IL

## Abstract

An UHPLC-ESI/MS/MS method was developed for simultaneous quantification and confirmation of meclofenamic acid (MFA) in equine plasma. MFA was recovered from equine plasma by strong anion exchange solid phase extraction (SAX-SPE). The concentration of MFA was determined between 5 and 250 ng/mL by internal standard calibration using flunixin-d<sub>3</sub> as the internal standard. The limits of detection (LOD) were determined to be 1 ng/mL. The measurement precision and accuracy were 4.7% and 92.9%, respectively. The method is not interfered by any other NSAIDs that are regulated by the United State Equestrian Federation (USEF). The method can assist USEF to control doping in horse racing.

## Introduction

MFA is a non-steroidal anti-inflammatory drug (NSAID). It is used to treat acute and chronic inflammation of the musculoskeletal system, including soft-tissue injury, bone and joint pathology, and laminitis. It is thought to be particularly useful in chronic problems of the hoof, including navicular syndrome and laminitis [1]. In racehorses, however, MFA has the capacity to affect racing performance, thus its usage is subject to doping control regulations. In order to discriminate between doping and therapeutic use of MFA, a threshold of 2.5 µg/mL in equine plasma has been adapted in samples collected immediately after the race by the United States Equestrian Federation (USEF).

Post-race analysis of MFA in equine plasma usually begins with a presumptive test using an enzyme-linked immunosorbent assay (ELISA). Positive samples are further submitted for quantification and confirmation of MFA by chromatographic based methods. For the quantification of MFA in equine plasma, we have developed a liquid chromatography ultraviolet detection (LC-UV) method after MFA was recovered from equine plasma by SAX-SPE. However, LC-UV lacks defensible confirmation of analytes, so it will not withstand legal challenges in court. In contrast, liquid chromatography mass spectrometry (LC-MS) can extract molecular fingerprints from very low concentrations of analytes, so it can provide legally defensible confirmation and will withstand legal challenges in court [2]. In this study, we would like to develop a LC-MS method for the simultaneous quantification and confirmation of MFA in equine plasma.

## Experimental

### SAX-SPE procedure

- Sample pretreatment:** Dilute 200 µL equine plasma and 200 µL internal standard solution with 1.6 mL of 0.5% (V/V) ammonia in water.
- Column conditioning:** Use gravity flow to apply 1 mL methanol, followed by 1 ml HPLC water through the column.
- Sample loading:** Use gravity flow to load the sample.
- Column wash:** Use gravity flow to apply 2×500 µL water, followed by 2×500 µL methanol. Finally, use a gentle vacuum to dry the column.
- Analyte elution:** Use gravity flow to apply 2×500 µL 5% (V/V) formic acid in methanol.
- Eluate drying:** Use gentle nitrogen stream to dry the eluates.
- Sample resuspension:** Add 200 µL mobile phase.

Table 1. Agilent 1260 Infinity II LC conditions

Parameter	Value
Column	Agilent Zorbax Eclipse C18 50 mm × 2.1 mm, 1.8 µm
Column temperature	40 °C
Injection volume	10 µL
Mobile phase	A: Water/acetonitrile 95/5 + 0.1 % formic acid B: Methanol/acetonitrile 95/5
Flow rate	0.300 mL/min
Gradient program	
0.0 minute	40% B
2.0 minute	90% B
6.0 minute	90% B
6.5 minute	40% B
Stop time	6.5 minute
Post time	5.5 minutes

Table 2. Agilent 6545 Q-TOF MS and MS/MS parameters

Parameter	Value
System tune	Standard 3200 m/z; 2 GHz Extended dynamic range; high resolution slicer mode
Transmission tune	50–750 m/z; 2 GHz Extended dynamic range; high sensitivity slicer mode
Mass calibration	50–750 m/z; 2 GHz Extended dynamic range; high sensitivity slicer mode
Ion source	Dual AJS ESI
MS acquisition mass range	100–1000 m/z
MS acquisition rate	10 spectra/s
MS/MS acquisition mass range	50–350 m/z
MS/MS acquisition rate	10 spectra/s
Drying gas temperature	325 °C
Drying gas flow	10 L/min
Nebulizer pressure	20 psi
Sheath gas temperature	400 °C
Sheath gas flow	12 L/min
Ionization mode	Positive
Capillary voltage	3500 V
Nozzle voltage	800 V
Fragmentor	120 V
Skimmer	45 V
Oct1 RF Vpp	750 V
MS reference mass ions	121.0509, 922.0098

Table 3. Agilent 6545 Q-TOF MS and MS/MS parameters

Parameter	Precursor (m/z)	RT (min)	Δ <sub>RT</sub> (min)	Isolation width	CE	Quantifier ion	Quanlifier ion
MFA	296.024	5.62	0.5	~1.3 m/z	20	243.045	278.009
Flunixin-d <sub>3</sub>	300.103	4.88	0.5	~1.3 m/z	25	282.293	264.050

## Results

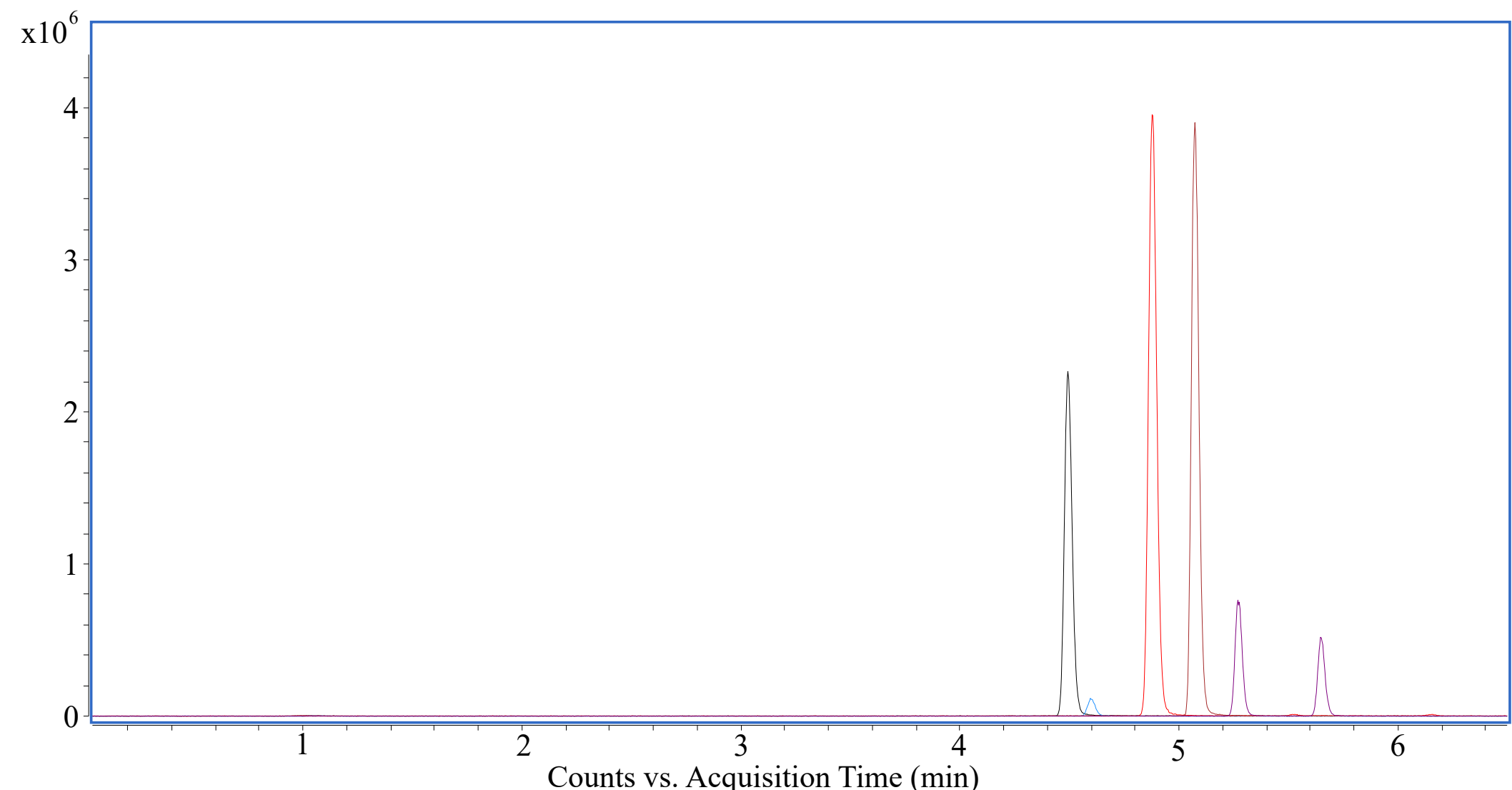


Figure 1. UHPLC-ESI/MS baseline separation of MFA from other NSAIDs regulated by USEF. The peaks were identified by their retention time as listed in Table 4. Each peak was extracted using the [M+H]<sup>+</sup> ion of the analyte (Table 4) with ±20 ppm mass window

Table 4. UHPLC-ESI/MS interference study: retention time, theoretically calculated *m/z*, and experimentally measured *m/z* of NSAIDs regulated by USEF

NSAID	Retention time (minutes)	[M + H] <sup>+</sup> (Calculated)	[M + H] <sup>+</sup> (Measured)
Ketoprofen	4.49	255.1015	255.1016
Naproxen	4.60	231.1016	231.1016
Flunixin	4.88	297.0841	297.0845
Phenylbutazone	5.07	309.1595	309.1598
Diclofenac	5.27	296.0240	296.0240
MFA	5.65	296.0238	296.0240

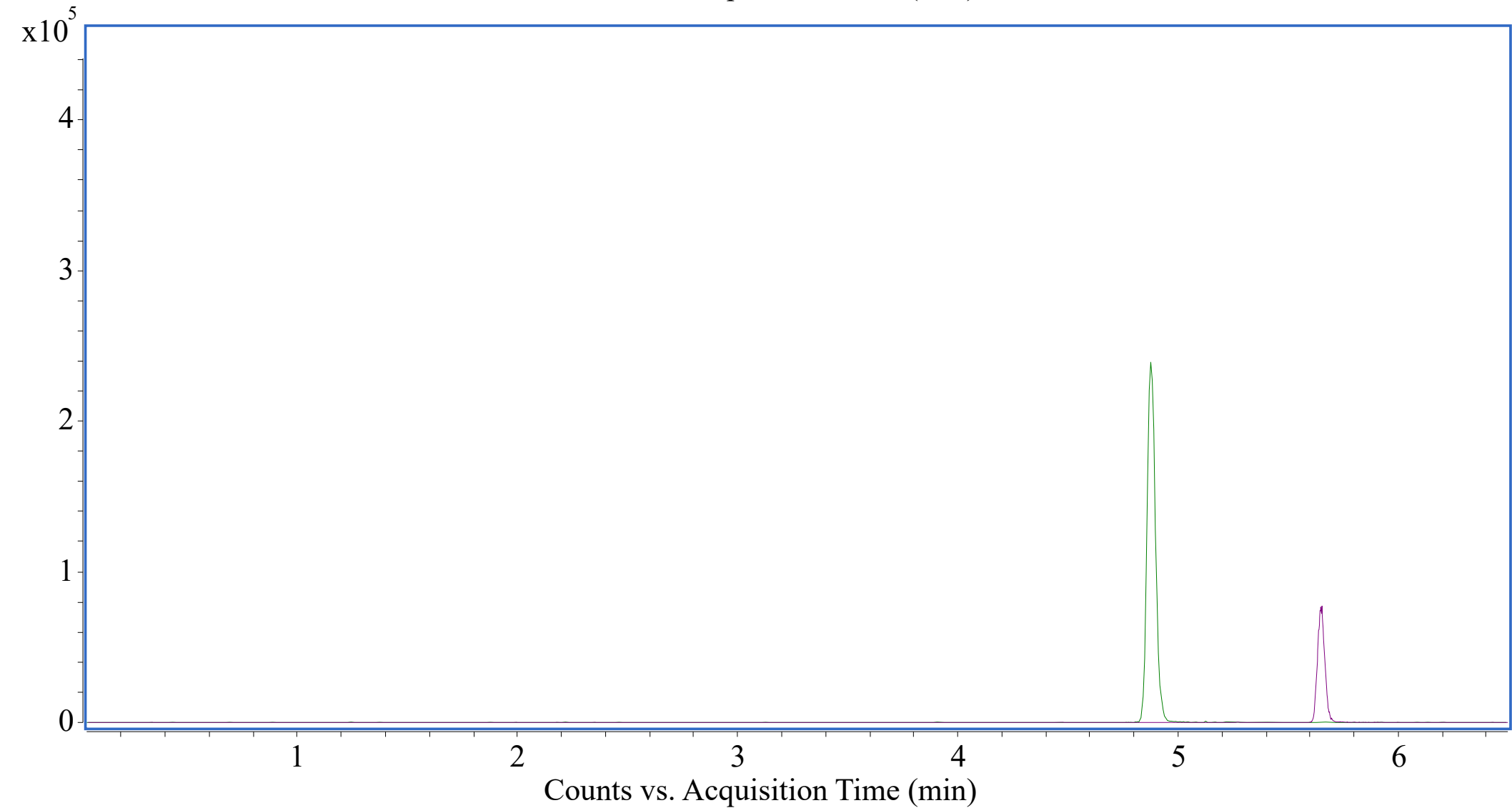
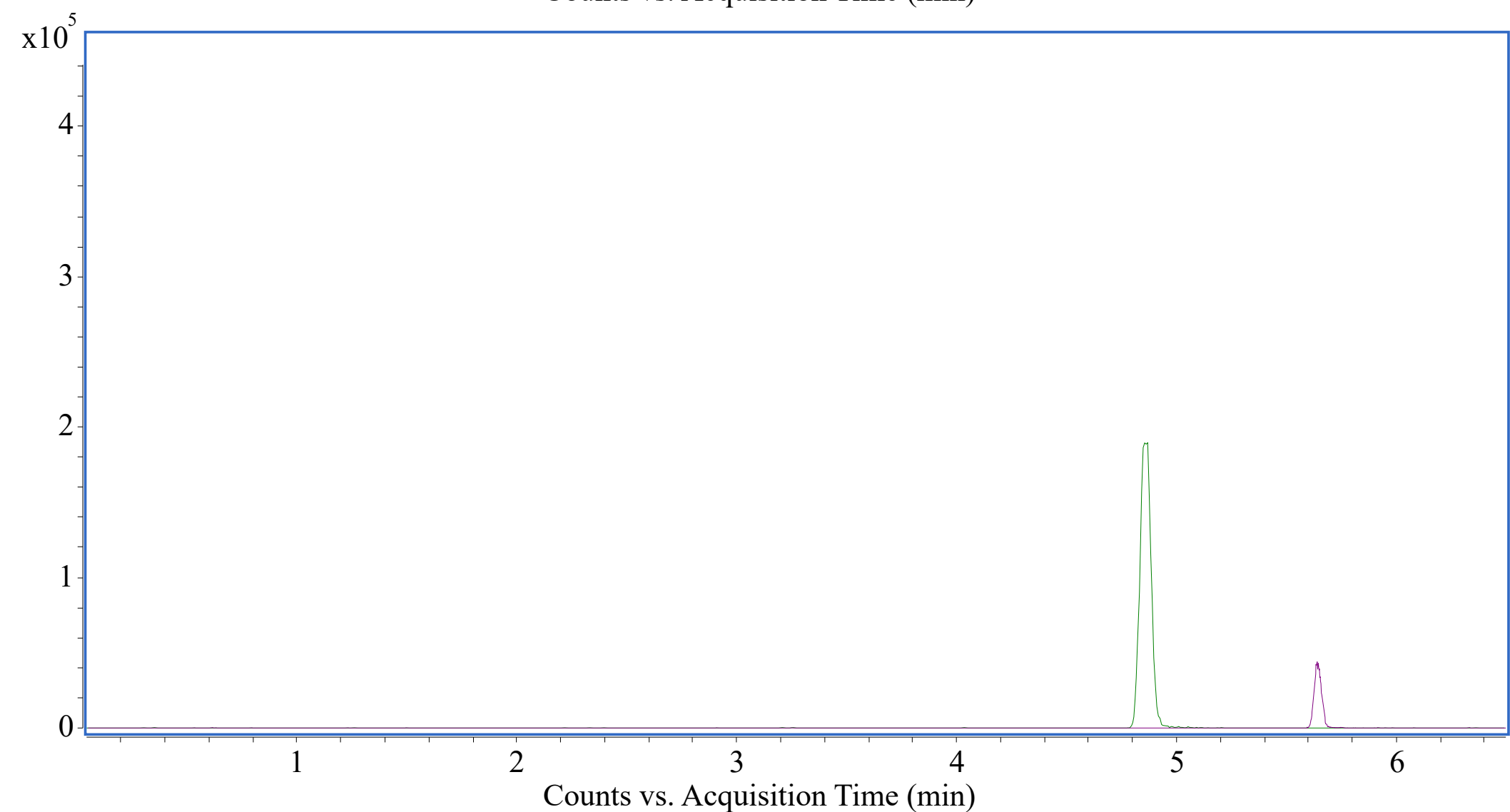
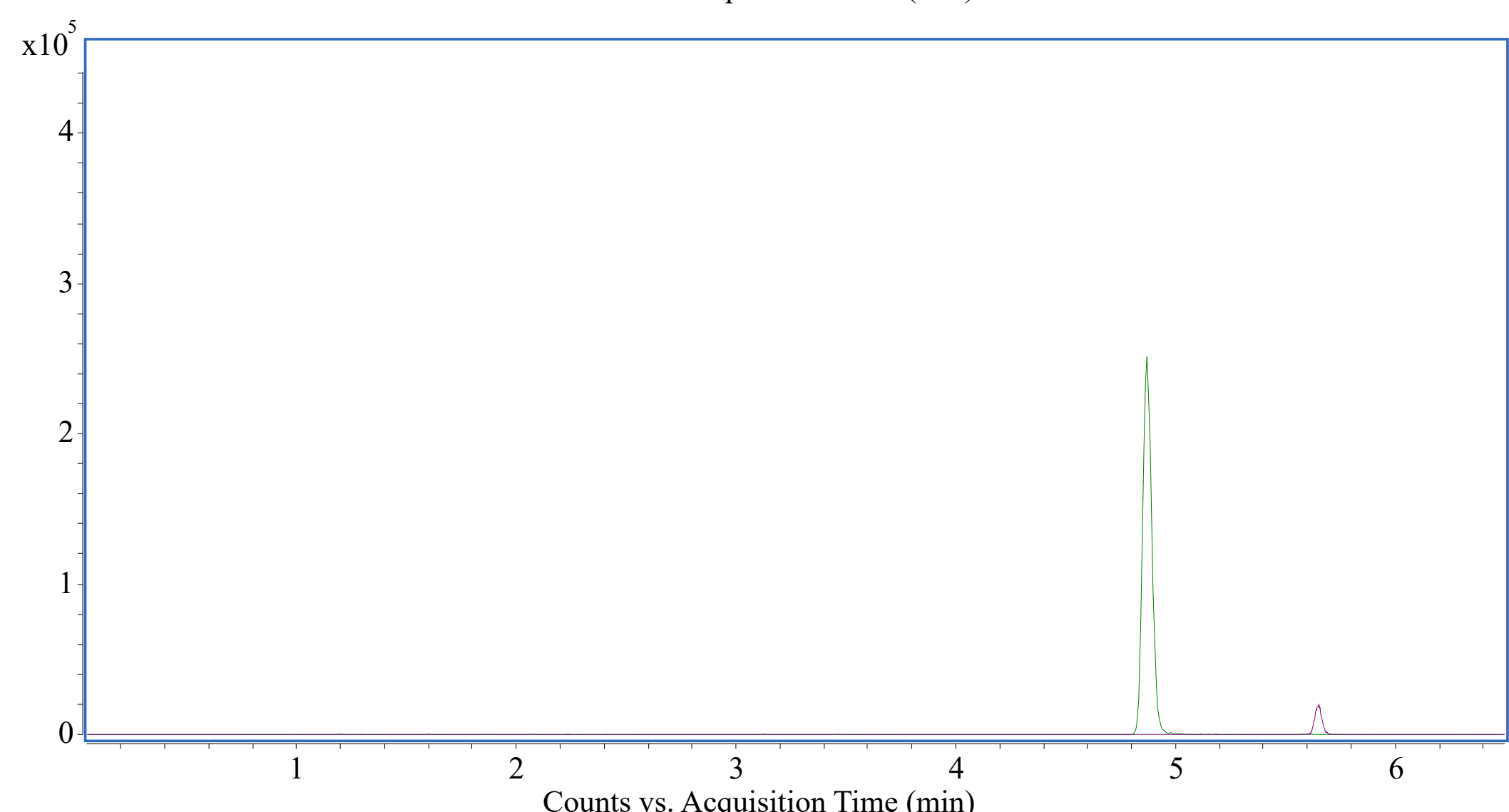
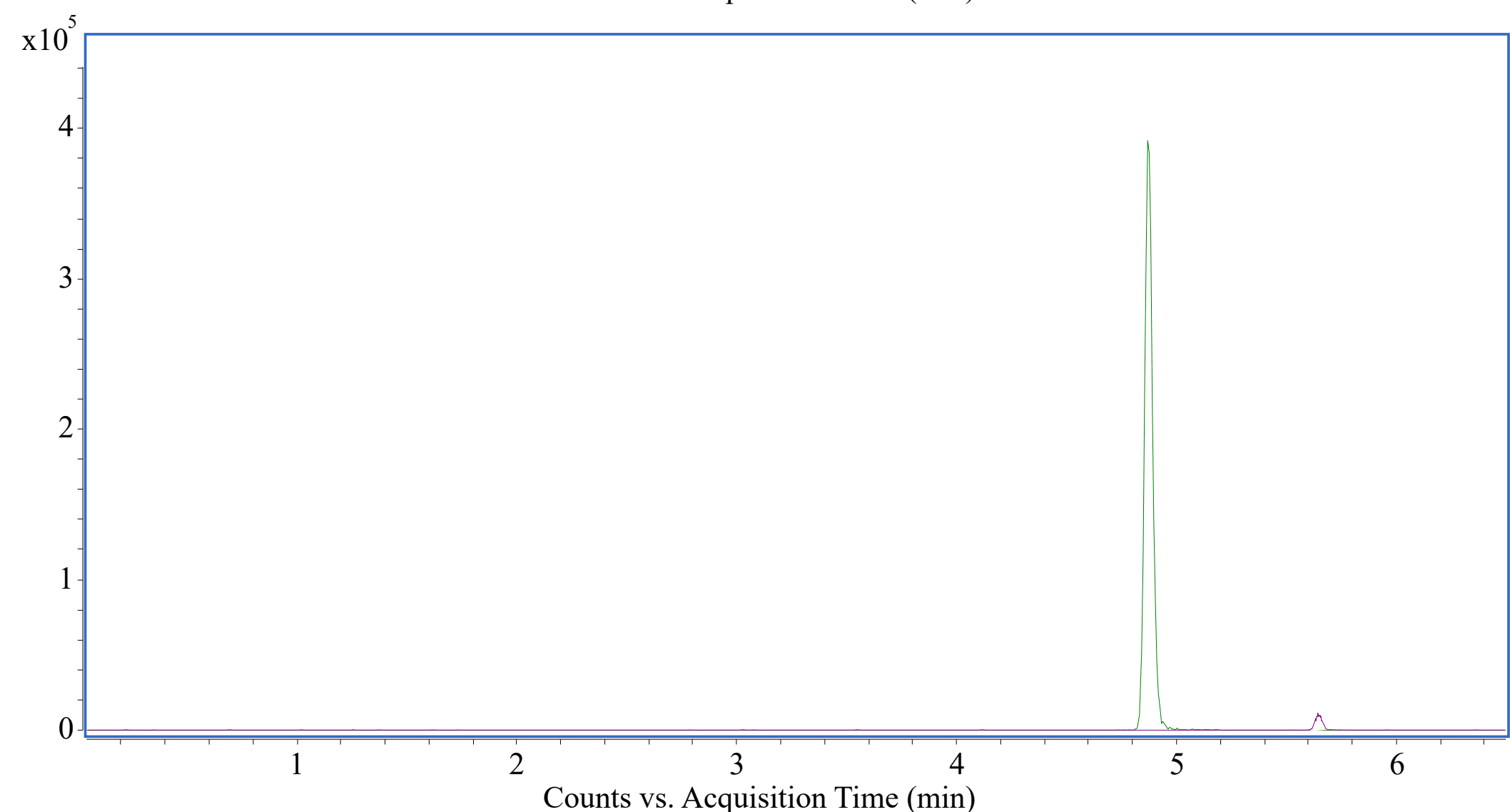
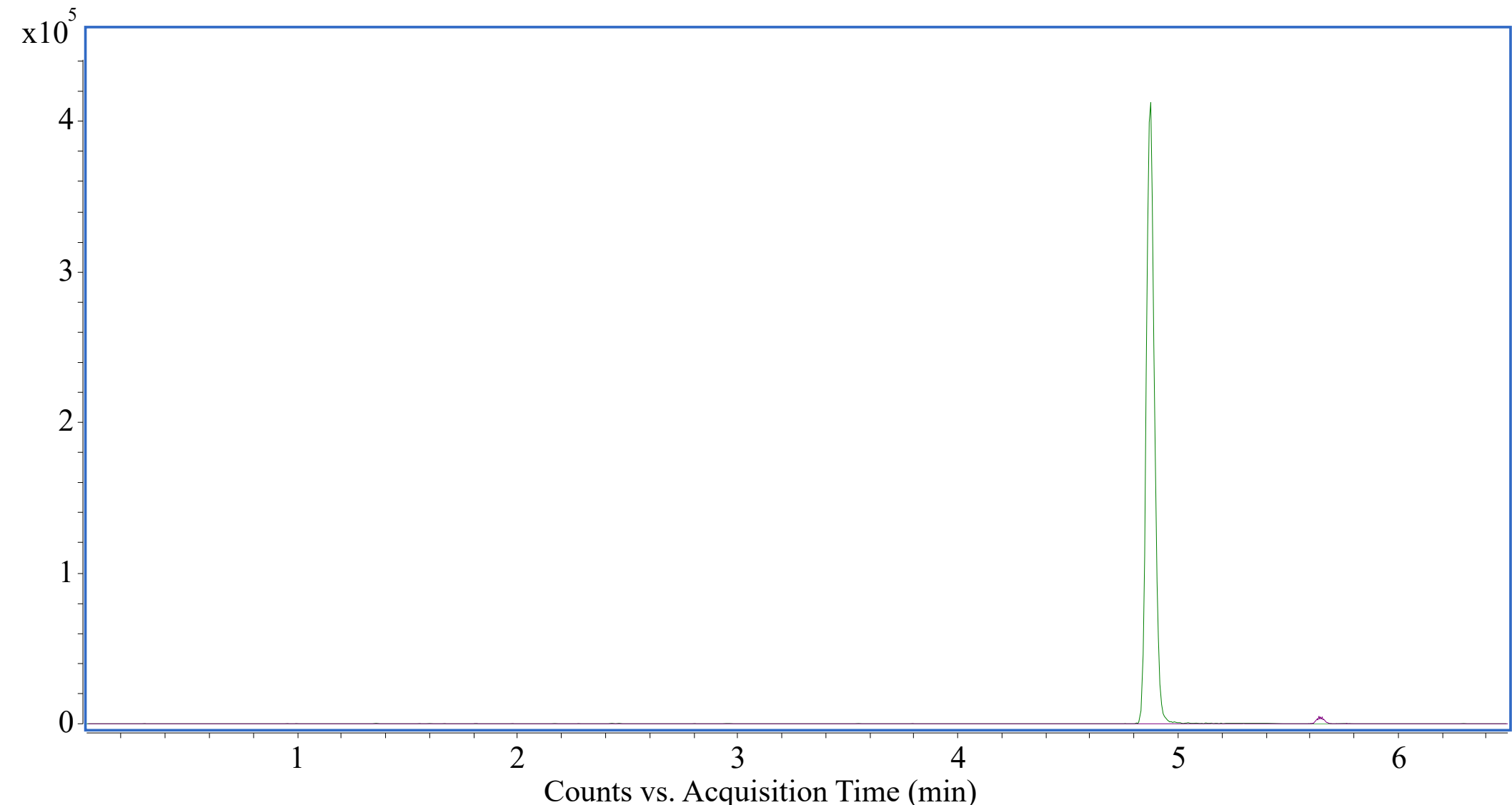


Figure 2. UHPLC-ESI/MS/MS chromatogram of 10, 25, 50, 125 and 250 (from top to bottom) ng/mL MFA (purple) and 25 ng/mL flunixin-d<sub>3</sub> (green).

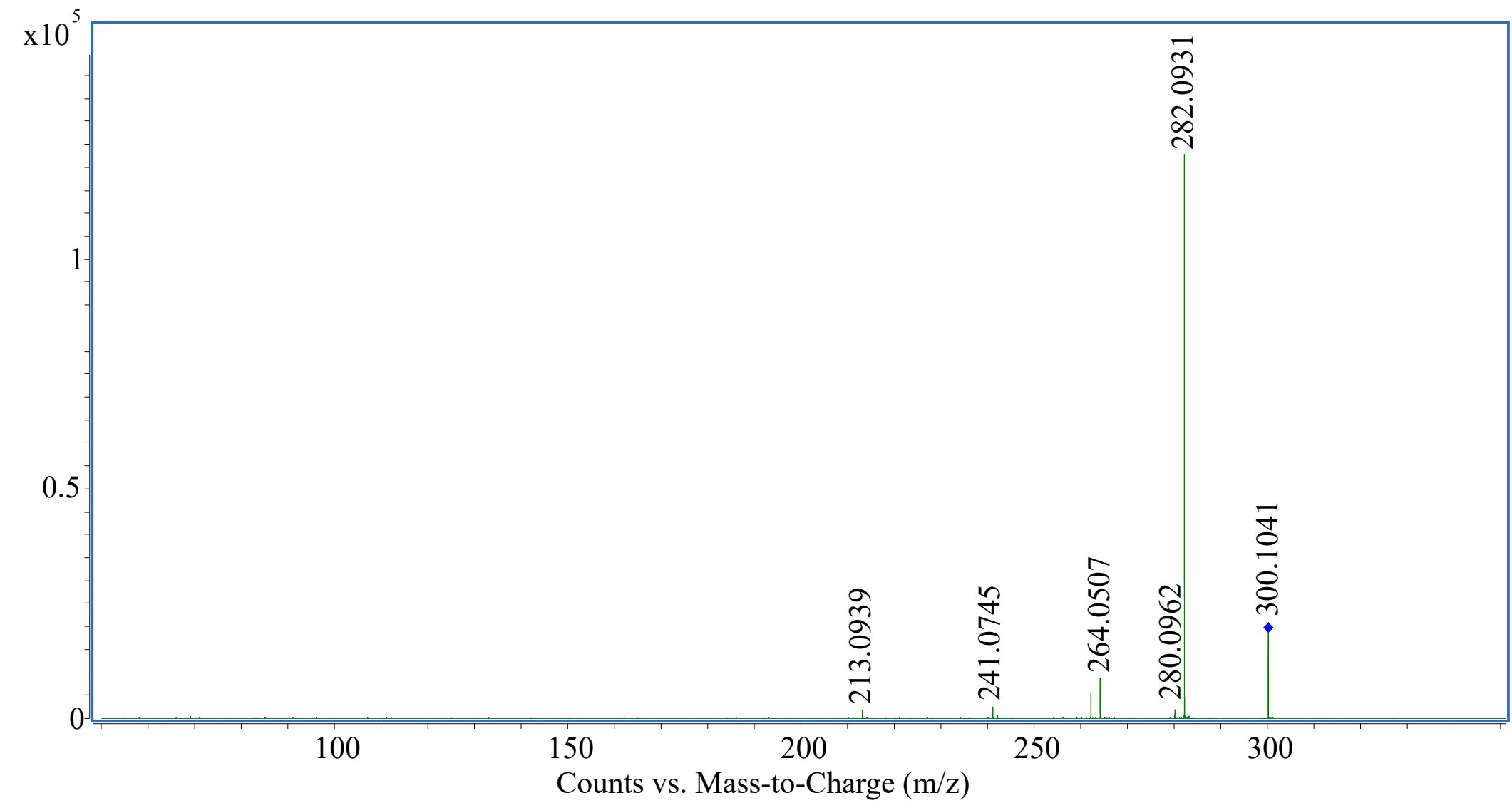
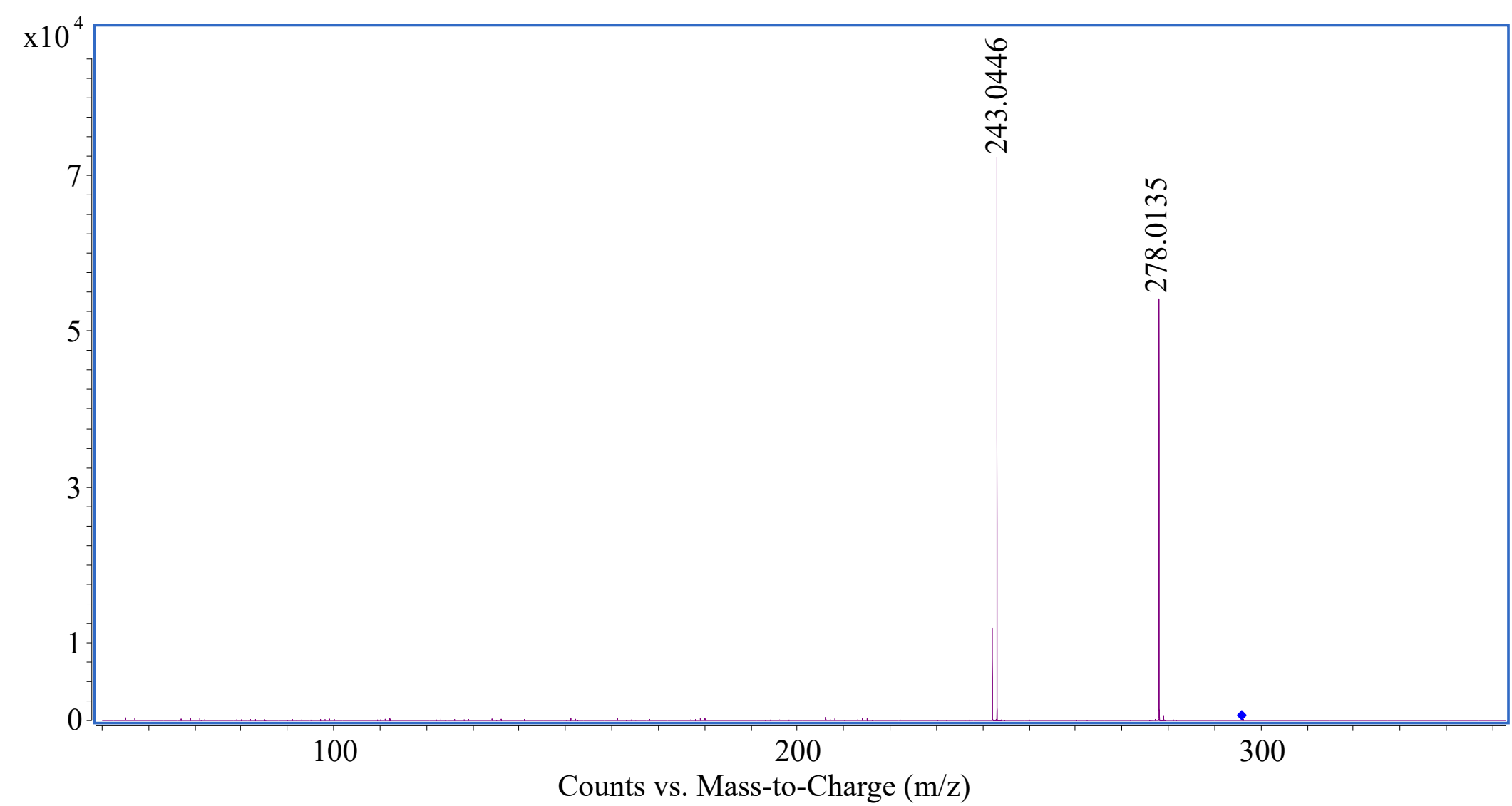


Figure 3. UHPLC-ESI/MS/MS spectra of MFA (purple) and flunixin-d<sub>3</sub> (green).

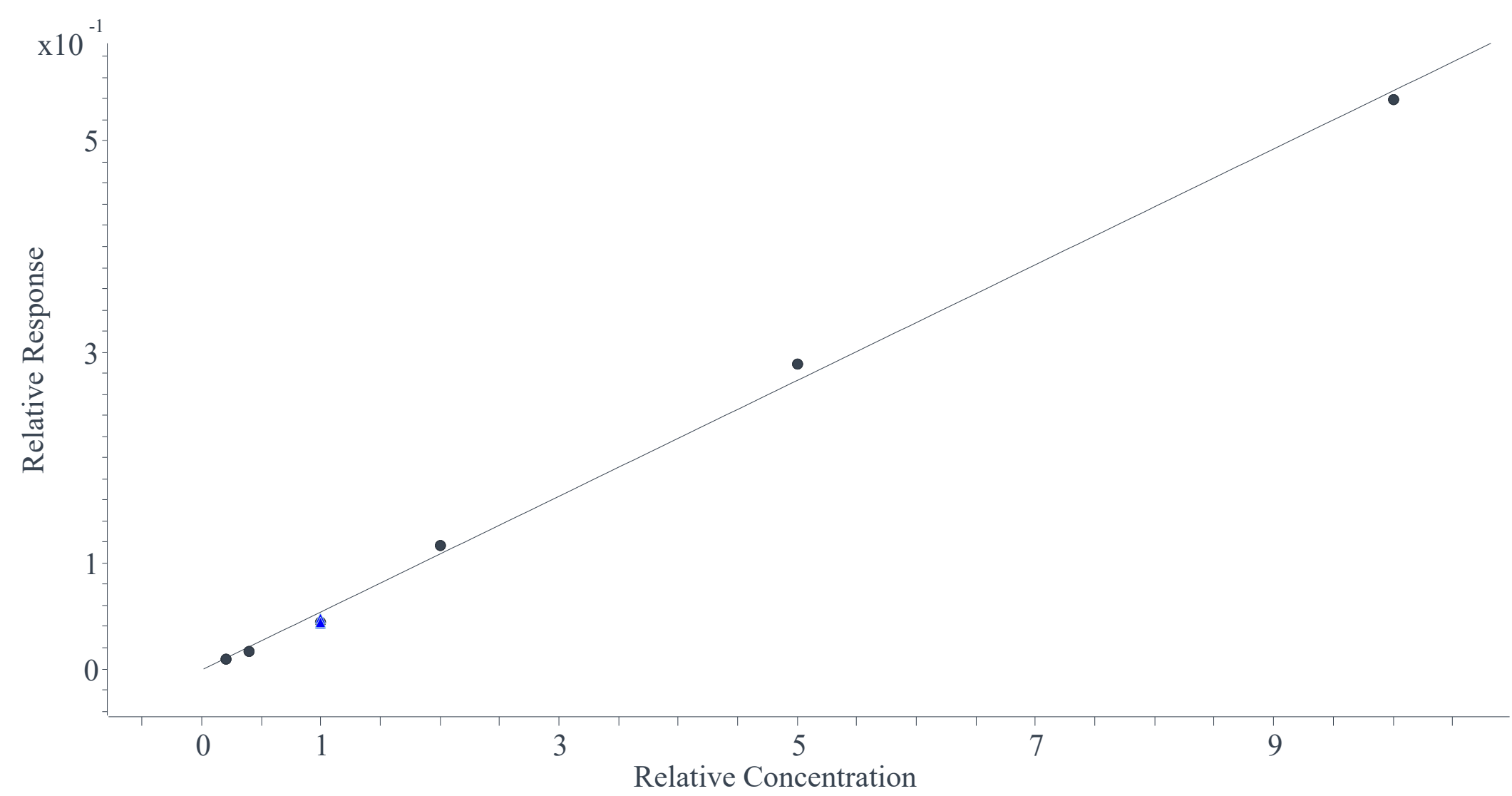


Figure 4. MFA internal standard calibration curve

Table 5. Precision and accuracy

QC	Spiked Conc.(ng/mL)	Measured Conc. (µg/mL)	% Accuracy
1	25.0	22.0	88.0
2	25.0	24.1	96.2
3	25.0	23.6	94.5
Average		23.2	92.9
% RSD		4.7	4.7

## Conclusions

- An UHPLC-ESI/MS/MS method for simultaneous quantification and confirmation of MFA in equine plasma has been developed.
- The limits of detection (LOD) was 1 ng/mL.
- The precision and accuracy were 4.8% and 101.9%, respectively.
- The method is not interfered by any other NSAIDs that are regulated by the United State Equestrian Federation (USEF).
- The method can assist USEF to control doping in horse racing.

## References

- Knych, H.K. Vet. Clin. N. Am.-Equine Pract. 33, 1-15 (2017).
- B. Heffron et al. J. Anal. Toxicol., 37 (2013) 600-604.

## Acknowledgement

The authors appreciate the financial support by the National Science Foundation (NSF) (CHE-1827209) and the Illinois State Academy of Science (ISAS).