

SLIDE: Saliva-based SARS-CoV-2 Self-Testing with RT-LAMP In a Mobile Device

Zifan Tang¹, Jiarui Cui¹, Aneesh Kshirsagar¹, Tianyi Liu¹, Michele Yon², Suresh V. Kuchipudi^{2, 3}, Weihua Guan^{1, 4*}

¹ Department of Electrical Engineering, Pennsylvania State University, University Park 16802, USA.

² Animal Diagnostic Laboratory, Pennsylvania State University, University Park, Pennsylvania 16802, United States

³ Center for Infectious Disease Dynamic, Pennsylvania State University, University Park, Pennsylvania 16802, United States

⁴ Department of Biomedical Engineering, Pennsylvania State University, University Park 16802, USA.

* Correspondence should be addressed to wzg111@psu.edu

Abstract

Regular, accurate, rapid, and inexpensive self-testing for SARS-CoV-2 is urgently needed to quell pandemic propagation. Existing at-home nucleic acid testing (NAT) test has high sensitivity and specificity, but it requires users to mail the sample to the central lab, which often takes 3-5 days to obtain the results. On the other hand, rapid antigen tests for SARS-COV-2 antigen provide a fast sample to answer the test (15min). However, the sensitivity of antigen tests is 30% to 40% lower than nucleic acid testing, which could miss a significant portion of infected patients. Here, we developed a fully integrated SARS-CoV-2 reverse transcription loop-mediated isothermal amplification (RT-LAMP) device using a self-collected saliva sample. This platform can automatically handle the complexity and can perform the functions including (1) virus particles thermal lysis preparation, (2) sample dispensing, (3) target sequence RT-LAMP amplification, (4) real-time detection, and (5) result report and communication. With a turnaround time of less than 45 minutes, our device achieved the limit of detection (LoD) of 5 copies/ μ l of the saliva sample, which is comparable with the LoD (6 copies/ μ l) using FDA-proved qRT-PCR assays with the same heat-lysis saliva sample preparation method. With clinical samples, our platform showed a good agreement with the results from the gold standard RT-PCR method. These results show our platform can perform self-administrated SARS-Cov-2 nucleic acid testing by laypersons with noninvasive saliva samples. We believe our platform self-testing platform will have an ongoing benefit for COVID-19 control and fighting future pandemics.

Keywords: SARS-CoV-2, POCT, Self-testing, RT-LAMP, Saliva

Coronavirus disease 2019 (COVID-19) became a worldwide pandemic in early 2020 ¹, and it was rapidly announced as a public health emergency of international concern by the World Health Organization (WHO) ^{2,3}. As of March 2022, there are more than 400 million confirmed cases and 6 million deaths of SARS-CoV-2 reported globally ³. Due to the fast mutation nature of the RNA virus and so many asymptomatic cases, all countries still face an unmet need to achieve a rapid, sensitive and reliable way to tackle the global and urgent problem. So far, nucleic acid amplification test (NAAT), such as RT-PCR, is the gold-standard technique due to its high sensitivity and specificity ⁴⁻⁷. However, laboratory-based NAAT requires highly trained personnel, dedicated facilities, and instrumentations, which typically require 3-5 days to get the result. Moreover, taking the on-site test requires people to stay with other potential patients, increasing the exposure risk. To alleviate these bottlenecks, the COVID-19 home test has become a practical option. Two different COVID-19 home tests are available: at-home PCR test ⁸ and antigen rapid test (Ag-RDT) ^{9,10}. So far, the FDA has issued EUA COVID-19 home tests developed by LabCorp, EverlyWell, Quest Diagnostics, PrivaPath Diagnostics, and Clinical Reference Laboratory. The user can self-collect the sample and shipped overnight to the company laboratory for SARS-CoV-2 viral RNA detection by PCR. Results are usually provided to test subjects within 3 -5 days ⁸. Even though it decreased the exposure risk, but longer time to obtain the results will increase the virus spread and delay the treatment. Ag-RDT tests are fast and cheap. It identifies active infection by detecting SARS-CoV-2 viral proteins. From sample collection to result takes 15–20 min using a portable device ^{11,12}. But the sensitivity of antigen tests is typically 30% to 40% lower than the nucleic acid testing ^{10,13,14}. Especially after the acute phase, when the viral load decreases, Ag-RDT might lead to high rates of false negatives, which could miss a significant portion of infected patients ⁹.

To overcome the drawbacks of the long sample-to-result time of conventional NAAT and the less sensitive rapid antigen test, developing a home-used sample-in answer-out NAAT analyzer for rapid and accurate COVID-19 detection becomes extremely necessary. Molecular diagnostics typically has five essential steps: (1) cells or virus particles lysis and DNA or RNA extraction, (2) sample partition, (3) target sequence amplification, (4) real-time detection by optical or other types of sensing mechanism ^{15,16}, and (5) data processing and result report. Integrating all those functions in a single device is critical to achieving self-testing and speeding up the process.

Since August 2020, saliva has become an alternative sample type for SARS-CoV-2 detection

¹⁷⁻¹⁹. This easy, noninvasive method largely increases the accessibility of self-testing ^{5, 20-24}. Several studies have demonstrated that saliva has comparable performance with nasopharyngeal samples ²⁵. Moreover, this saliva sample preparation has been further simplified by the Yale school of public health researcher. They found that 5 minutes of heat-inactivation of the saliva sample without any additional reagents can achieve a low limit of detection (6 copies/ μ L) using FDA-approved RT-PCR assays ²¹.

Recently, isothermal amplification techniques have been widely used for the point of care setting, for example, reverse transcription loop-mediated isothermal amplification (RT-LAMP) ^{5, 26-42}. The RT-LAMP process is similar to conventional PCR tests, but the reaction can be performed without commercial thermocyclers. While maintaining specificity and sensitivity comparable to that of the PCR tests, RT-LAMP shows better tolerance for the impurities and a faster time to result. These unique features make RT-LAMP assays quicker, easier to use, and more cost-effective than RT-PCR assays, making them more suitable for POC diagnostics.

In this work, we developed a fully integrated SARS-CoV-2 NAT device using a self-collected saliva sample. This Saliva-based SARS-Cov-2 Self-Testing with RT-LAMP In a Mobile Device (SLIDE) platform consists of a ready-to-use reagents cartridge, an easy-to-use smartphone interface, and an ultra-compact analyzer. It automatically handled the complexity of heat-inactivated sample preparation, sample dispensing, real-time RT-LAMP reaction and detection, and data processing. With a turnaround time of less than 45 minutes, we achieved a limit of detection (LoD) of 5 virion/ μ L of a saliva sample. With clinical samples, our platform showed a good agreement with the results from the gold standard RT-PCR method. We believe our platform self-testing platform will have an ongoing benefit for COVID-19 control and fighting future pandemics.

Results and Discussion

Overall Design and Module Validations

Overall Design. The overall design of the SLIDE analyzer is shown in **Figure 1a**. It consists of five seamlessly integrated modules controlled by a microcontroller unit (MCU): an optical module for excitation and detection, two thermal modules, a piezo micro pump module, a power

module, and a connectivity module. **Supplementary Figure S1** illustrates the overall block diagram design of the device. The whole system is designed in SolidWorks and prototyped with in-house 3D printing. **Figure 1b** shows a photograph of the assembled SLIDE analyzer and the smartphone interface.

Optical Module. The optical module consists of three independent excitation and detection units. Each unit has a LED excitation source ($\lambda=470$ nm) and a CMOS color sensor for real-time fluorescence monitoring. The excitation and the detection were arranged to be perpendicular to each other to minimize the excitation interference on the fluorescence signal (**Figure 1c**). To characterize the quantification ability of the optical module, we tested different calcein concentrations from 0 to 25 μM and measured the fluorescence intensity for 10 minutes. **Figure 1d** showed the mean and standard deviation of the relative fluorescence unit (RFU) as a function of the calcein concentration. A linear fit with $R^2=0.98$ confirmed the quantitative capability of the optical module.

Thermal Module. We designed two independent thermal modules. One is for heat-inactivating the saliva and performing the thermal lysis at 95°C. The other is for controlling the temperature of the RT-LAMP reaction at 65°C. Both modules used a customized aluminum heating block with power resistors attached. The temperature was controlled through a feedback measurement of a thermistor embedded in the heating block. Since the temperature was obtained from the heating block rather than the analyte solution on the cartridge, we characterized the temperature profile difference between these two. As shown in **Figure 1e**, the heating block reached 95 °C after 2 mins of operation, while the saliva in the cartridge took 5 mins. This delay is due to the non-ideal thermal coupling and the different specific heat capacity between the heating block and the cartridge. Nevertheless, the saliva can be sufficiently lysed at 95°C within 10 mins from sample collection. For the heating module controlling the RT-LAMP reaction, we observed that the mean and the standard deviation of the temperature in the master mix solution is 64°C and 0.38°C, respectively (**Figure 1f**).

Micro pump Module. The sample dispensing and mixing is accomplished on the cartridge using a micro piezo pump. It is connected to the microfluidic cartridge using a Tygon tube and a Luer-lock adaptor. The volumetric rate of the micro piezo pump is controlled by the frequency and the driving voltage. To characterize the micro pump, we tested the volumetric rate at different frequencies. As expected, the volumetric rate increased linearly with the operation frequency ($R^2=$

0.99, **Figure 1g**). This relationship provides us with the capability to modulate the liquid flow rate on the cartridge through programming the operation frequency.

Power Module. A rechargeable 1300 mAh Lithium polymer battery (14.43 Wh) was used to power our analyzer. To estimate the power consumption for each run, we used a power meter to characterize the voltage, current, and power during a complete cycle of the test. **Figure 1h** shows a complete-time trace. As shown, heating is the most power-hungry process during the operation. Before reaching the target temperature, the heaters continuously work at a high current (1.7 A for 95 °C and 2.2 A for 65 °C). After reaching the target temperature, the heater starts to change states between on and off to maintain the temperature. The total energy consumed is 3.02 Wh in each 45 minutes test, meaning we can perform at least four tests before recharging.

Connectivity Module. A smartphone app was developed to assist the user in conducting the test. The flow chart of the app process is shown in **Supplementary Figure S2**. The SLIDE analyzer and the smartphone communicated through Bluetooth LE protocol. The App could provide test instructions, acquire data, and make positive and negative calls to interpret the test results. The App could also save the test results into a spreadsheet, save them on the local smartphone, or upload them to cloud-based storage (Google Drive). **Supplementary Figure S3** shows the representative screenshots of the developed App.

Automated Salvia Processing on The Cartridge

To facilitate the raw salvia processing, we developed a disposable cartridge with the SLIDE analyzer. The cartridge was fabricated in polymethyl methacrylate (PMMA). It consists of three laminated layers: top layer, middle microchannel layer, and bottom layer (**Figure 2a**). The overall layout of the assembled cartridge shows in **Figure 2b**. It includes a heat lysis chamber (250 μ l), three independent dispensing (10 μ l) and reaction chambers (60 μ l), and a waste chamber (300 μ l). First, the collected raw saliva sample was heat-inactivated and lysed at 95°C for 5 minutes. The resulting lysates were transferred to the dispensing and reaction chambers through the microchannel. The excessive analyte sample was stored in the wastes chamber with a venting hole to the atmospheric pressure.

Figure 2c illustrates the detailed design of a single unit of dispensing and reaction chambers. Since the laser-processed PMMA side walls are hydrophilic ⁴³, a side pocket structure can easily

trap 10 μ l of the samples without bubbles. We found that the 5.3 ml/min flow rate could help ensure the reliable trapping process. During the trapping process, the paraffin wax valve 1 was in the solid phase such that the trapping volume was fixed (**Figure 2d**). The average trapping volume is 10.25 ± 0.27 μ l. The difference between the three chambers was less than 2.5% (**Supplementary Figure S4**).

After excessive samples were pushed into the waste chamber and each unit metered 10 μ l of the heat-processed saliva, we increased the temperature to 65 °C to open the paraffin wax valves. When the wax valves 1 and 2 change from the solid to the liquid phase, the trapped saliva sample will start flowing into the reaction chamber by capillary force. To facilitate transferring all samples into the reaction chamber and thorough mixing with the RT-LAMP master mix, we applied 30 consecutive micro pump pressure pulses. Each pulse is programmed to be 100 ms in duration (**Figure 2e**). The paraffin wax valve 2 serves as a hydraulic resistor, which helps to balance the hydraulic resistance among three units. To avoid liquid overflowing, we intentionally designed a long S-shaped releasing channel with a venting hole at the end. In addition, a thin layer of wax on top of the RT-LAMP mix protects the master mix from evaporation. It also avoids external contamination by providing a barrier against amplicons from escaping. **Figure 2f** and **Supplementary Video S1** showed a representative example of automated saliva processing on the cartridge.

Saliva Test Workflow

The overall SLIDE workflow from the saliva sample to the molecular results is shown in **Figure 3a**. Four components are needed for a test: a disposable cartridge, a saliva collection aid (SCA), a portable analyzer, and an Android smartphone. With the help of the instructions on an interactive smartphone app (**Figure 3b**), one would self-collect saliva samples into a cartridge with the help of an SCA. While collecting the whole saliva through spitting or drooling is feasible, the saliva collection aid could increase participant compliance and avoid sample foaming²¹. After sufficient saliva (~120 μ l) was collected into the cartridge, the user should seal the cartridge with a screw cap. The sealed cartridge can then be connected to the piezo pump through a Luer-lock interface and be inserted into the analyzer. One then would need to turn on the analyzer for the smartphone to recognize and communicate through the Bluetooth connection. This process takes

less than 2 min hands-on time and is the only manual testing step.

Once the SLIDE analyzer receives the ‘start testing’ command from the smartphone app, the analyzer will automatically perform the required tasks on the cartridge. It includes saliva thermal lysis, sample metering and dispensing, RT-LAMP reaction and real-time detection, and data analysis and storage. Specifically, the analyzer begins the test by thermal lysis of the saliva sample at 95°C for 5 minutes. This step inactivates RNases and releases the virus from the saliva sample²¹. The resulting lysates were automatically transferred and dispensed into the reaction chamber with a preloaded RT-LAMP master. The whole sample preparation takes about 13 minutes. After dispensing the sample, the real-time RT-LAMP reaction starts at a constant temperature of ~64 °C⁴¹. The acquired fluorescence data are transmitted to the smartphone app every 5 seconds. The threshold to distinguish the positive from the negative was set at 50 RFU based on the NTC samples tested (**Supplementary Figure S5**). We classify a sample as positive only when two out of three reactions have a higher RFU than the threshold value in 30 minutes. The test results could be saved on the local device and uploaded to a cloud. The whole process is fully automated (**Supplementary Figure S6**) and takes about less than 45 min (~2 min hands-on time for sample collection, ~13 min for sample preparation and dispensing, and ~15-30 min for RT-LAMP reaction, data processing and result report) from the saliva collection to the result, with very minimal user intervention (**Supplementary Video S2**).

Performance Evaluation with Mock Saliva Sample

After validating all the subsystems and system integration, we went out to test the performance of the SLIDE. Here, we used our previously validated SARS-CoV-2 RT-LAMP primer set⁴¹ (**Supplementary Table S1**) against the highly conserved N region with a modified fluorescent concentration of SYTO9 (**Supplementary Table S2**). We formed mock SARS-CoV-2 positive samples by spiking the healthy saliva with different concentrations of heat-inactivated SARS-CoV-2 virus particles. The final viral concentration of the mock sample ranges from 1 to 10⁴ copies/μl. **Figure 4a** shows the real-time result. Note that each sample is aliquoted to three separate reactions on a single cartridge (**Figure 2b**). The sample is classified as positive in each test only if more than two out of three reactions have an RFU of more than a threshold. As shown, samples with a concentration above 5 copies/μl were successfully classified as positive, while 1 out of 3

samples at 1 copy/ μ l were classified as positive.

To estimate the LoD of the test, we examine the hit rates at different virus concentrations⁴⁴. The hit rate is the positive test over all the tests under the same concentration. As shown in **Figure 4b**, the hit rate started to roll off from 100% to 33% when the concentration decreased from 5 copies/ μ l to 1 copy/ μ l. We fitted the experimental hit rate data with a logistic curve. The LoD is determined to be about 5 copies/ μ l at the 98% confidence level. This LoD is comparable with the LoD (6 copies/ μ l) using FDA-proved qRT-PCR assays with the same heat-lysis saliva sample preparation method²¹. **Figure 4c** shows the threshold time in the SLIDE analyzer with different virus concentrations. The threshold time and the standard deviation among the times to positive generally increases as the virus particle concentration decreases, although the linearity is not as good as a RT-PCR test. The less ideal linearity is expected as RT-LAMP assay is not a quantitative assay.

To further evaluate our device, the same spike samples were tested using the benchtop PCR machine. We manually performed the sample thermal lysis in the heating block for 5 min at 95°C, then transferred 10 μ l of the processed sample using a pipette to the PCR tube with a preloaded RT-LAMP master mix. After mixing the reagents thoroughly, the reactions were performed using a benchtop PCR machine (**Supplementary Figure S7**). **Figure 4d** show a Pearson correlation of the threshold time between the SLIDE analyzer and the PCR instrument. A coefficient ($R=0.835$) indicates a good agreement between the automated SLIDE device and manual methods.

Clinical Saliva Sample Test

To best evaluate the performance of SLIDE, clinical samples are tested. Here, two archived clinical samples (one known positive and one known negative) were obtained through an approved institutional review board (IRB) of the Pennsylvania State University. All the samples were coded to remove information associated with patient identifiers. The RT-PCR assay performed the initial diagnosis as the reference method to benchmark our SLIDE. The experiment follows the protocols shown in **Figure 3** and **Supplementary Video S2**. The resulting raw amplification curves are shown in **Figure 4e**. In 30 minutes of the amplification process, all three reactions in the positive test showed sharp RFU increases and stabilized at the RFU value at least three times above the threshold. All reactions in the negative clinical sample showed no noticeable RFU changes. The

positive and negative samples determined by the SLIDE analyzer agree with that of the gold standard RT-PCR method. Evaluating more clinical samples would be needed to demonstrate the device robustness and reproducibility. A scaled-up test with more clinical samples is currently under another IRB approval. We will evaluate and report the diagnostic sensitivity and specificity of the SLIDE device when these data are acquired in the future.

Conclusion

We demonstrated a fully integrated device for rapid (<45 min) self-testing of the SARS-CoV-2 virus from saliva samples. This fully portable device can detect the virus rapidly without needing an RNA extraction kit and pipetting steps. All other complexities are handled automatically by the SLIDE analyzer, including sample processing and dispensing, real-time RT-LAMP reaction and detection, and data processing and communication. Our automatic system shows an excellent agreement with the manual process using a benchtop PCR instrument. The limit of detection against the SARS-CoV-2 virus particle spiked in the saliva sample is 5 copies/ μ l. This LoD is comparable with the LoD (6copies/ μ l) using FDA-proved qRT-PCR assays with the same heat-lysis saliva sample preparation method ²¹. A pilot clinical saliva sample test with the SLIDE showed a good agreement with the gold standard RT-PCR method. These results show that it is feasible to perform self-administrated SARS-Cov-2 nucleic acid testing by laypersons with noninvasive saliva samples. To that end, we will need to further address the outstanding issues of reagent lyophilization on the cartridge and scaled up clinical validation in future studies.

Materials and methods

SARS-CoV-2 samples

Heat-inactivated SARS-CoV-2 (ATCC VR-1986HK) virus particle was purchased from ATCC. The negative saliva samples were collected from healthy volunteers. The mock samples were prepared by spiking the heat-inactivated SARS-CoV-2 virus particles into the healthy saliva sample. The pre-identified clinical saliva samples were approved by the institutional review board

(IRB). These clinical saliva samples were initially tested with the F.D.A. EUA-Authorized OPTI RT-PCR COVID-19 Direct assay (OPTI Medical Systems, GA, USA). The collected saliva samples were frozen at -80°C before use. All the clinical experiments were performed in the Animal Diagnostic Laboratory (BSL 3) at Penn State (University Park) by a protocol approved by the Institutional Biosafety Committee.

RT-LAMP reaction mix

The total volume (40 μ l) of the RT-LAMP assays contains a 30 μ l master mix and 10 μ l saliva sample. The master mix includes isothermal buffer, PCR grade H₂O, MgSO₄ (7 mM), Sty-9 green (10 μ M), deoxyribonucleotide triphosphates (dNTPs, 1.4 mM), Bst 2.0 DNA polymerase (0.4 U/ μ l), Warmstart reverse transcriptase (0.3 U/ μ l), primer sets (0.2 mM F3 and B3c, 1.6 mM FIP and BIP, 0.8 mM LF and LB, see **Supplementary Table S1** for primer design). (**Supplementary Table S2** summarizes the RT-LAMP recipe).

Instrumentation

Figure 1a shows a photo of a SLIDE analyzer. The SLIDE analyzer comprises 3D printed structural parts, a CNC machined aluminum heating block, a micro pump, electronics such as an Arduino Nano (MCU), excitation LEDs, color sensors for fluorescence detection and Bluetooth. The 3D housing was designed in Solidworks software and printed using MakerBot MethodX 3D printer (Brooklyn, NY) with MakerBot ABS (acrylonitrile butadiene styrene). The aluminum heating blocks were designed in Solidworks software and fabricated using a CNC machine. Two one-ohm power resistors are mounted (in series) on the aluminum heating using a thermally conductive adhesive paste for the 95°C heating block and 65°C heating block, respectively. Negative thermal feedback control was performed using N-channel power MOSFET (63J7707, Digi-Key) and an MC65F103A 10 k-ohm thermistor (Amphenol Thermometrics, St. Marys, PA) to maintain the desired temperature. PCBs were designed in AutoDesk Eagle CAD software and fabricated by O.S.H. Park L.L.C. (Lake Oswego, OR). The optical module PCB consists of three blue excitation LEDs (FD-5TB-1) purchased from Adafruit Industries (New York, NY) and three-color sensors (TCS 34725) purchased from DigiKey. The piezo pump and the driver were purchased from Bartels (Mikrotechnik, Germany). The Bluetooth (Adafruit Bluefruit LE SPI Friend) module was purchased from Adafruit Industries (New York, NY). The whole system was

powered by a 1300 mAh Lithium polymer battery (ZIPPY, USA). All materials of analyzer can be found in **Supplementary Table S3**.

APP development

Four steps are involved in this Android App development. First, the App interface guides users in providing their personal information. Only the name is required from users. The Global Positioning System (GPS) can automatically obtain the time and location information. Second, we set up Bluetooth communication. App interface scans and connects the Bluetooth LE around the analyzer. The communication protocol can be built using the Service UUID and Characteristic UUID of the Bluetooth LE, enabling the data communication function between these two devices. Once the user clicks the confirm button on the screen, the App will send a single bit to the analyzer to initiate the test. The third part is the real-time data transfer and plotting. We added two check bits at the beginning and the end of the string to ensure accuracy. After confirming the check bit of the received string from the analyzer, the string value will be split into three channels and plotted with different colors. Meanwhile, the split data in each channel is compared with the threshold value (RFU 50) to make the decision. If more than or equal to two channels have three successive data greater than the threshold, the test result will be identified as a positive. Otherwise, the App will continue to receive the string value from the analyzer. If no positive result has been determined after 30 minutes of the amplification, the test result will be negative. The APP will combine personal information, color sensor data in each channel, and test results into a spreadsheet. This file can be saved on the local device and uploaded to a Google drive. The flow chart of this APP development process is shown in **Supplementary Figure S2**. Selected screenshots of the APP are presented in **Supplementary Figure S3**.

Microfluidic reagent cartridge

The microfluidic cartridge was designed by AutoCAD and patterned using a CO₂ laser cutting machine (Universal Laser Systems, Scottsdale, AZ). All layers were aligned and laminated with an adhesive solvent (Weld-On). The assembled cartridge comprises a sample collection chamber (250 µl), three trapping chambers (10 µl each), three reaction chambers (60 µl each), three wax valves 1 (5 µl each), three wax valves 2 (5 µl each), and wastes chamber (300 µl). The sample collection tube was mounted using the Epoxy Adhesive (3M, Saint Paul, MN). All the assay and

wax valves are loaded onto the cartridge through the extruded inlet and sealed by the PCR plate seals (Bio-Rad, Hercules, CA). The Saliva Collection Aid is purchased from Salimetrics, LLC.

Data processing

To find the proper threshold, all the collected raw data is subtracted from the background signal acquired from the average of the first 10 data points and leveled at RFU 30. The threshold to classify an amplification curve as positive or negative was 50 RFU based on the negative sample (Supplementary Figure S5).

Associated Content

The Supporting Information is available. The analyzer system diagram, the flow chart of the android App development, mobile phone user interface, the trapping volume characterization, the representative measurement from negative control samples, the flow chart of the automatic workflow from sample to answer, real-time results using a PCR machine. Detailed descriptions of RT-LAMP primer set for N region, the recipe for RT-LAMP master mix, and the bill of material for the SLIDE analyzer.

Author contributions

WG conceived the concept and supervised the study. ZT, AK and TL developed the instrument and integration. ZT and JC developed the mobile APP. ZT developed and validated the RT-LAMP, cartridge design, and the device. ZT, SVK, MY performed the clinical samples study. ZT and WG co-wrote the manuscript, with discussion from all authors.

Conflicts of interest

A provisional patent related to the technology described herein is filed.

Acknowledgments

This work was supported by the National Institutes of Health (R61AI147419), National Science Foundation (1902503, 1912410, 2045169), and Coronavirus Research Seed Fund. Any opinions, findings, and conclusions or recommendations expressed in this work are those of the authors and do not necessarily reflect the views of the National Science Foundation and National Institutes of Health.

Figures and Captions

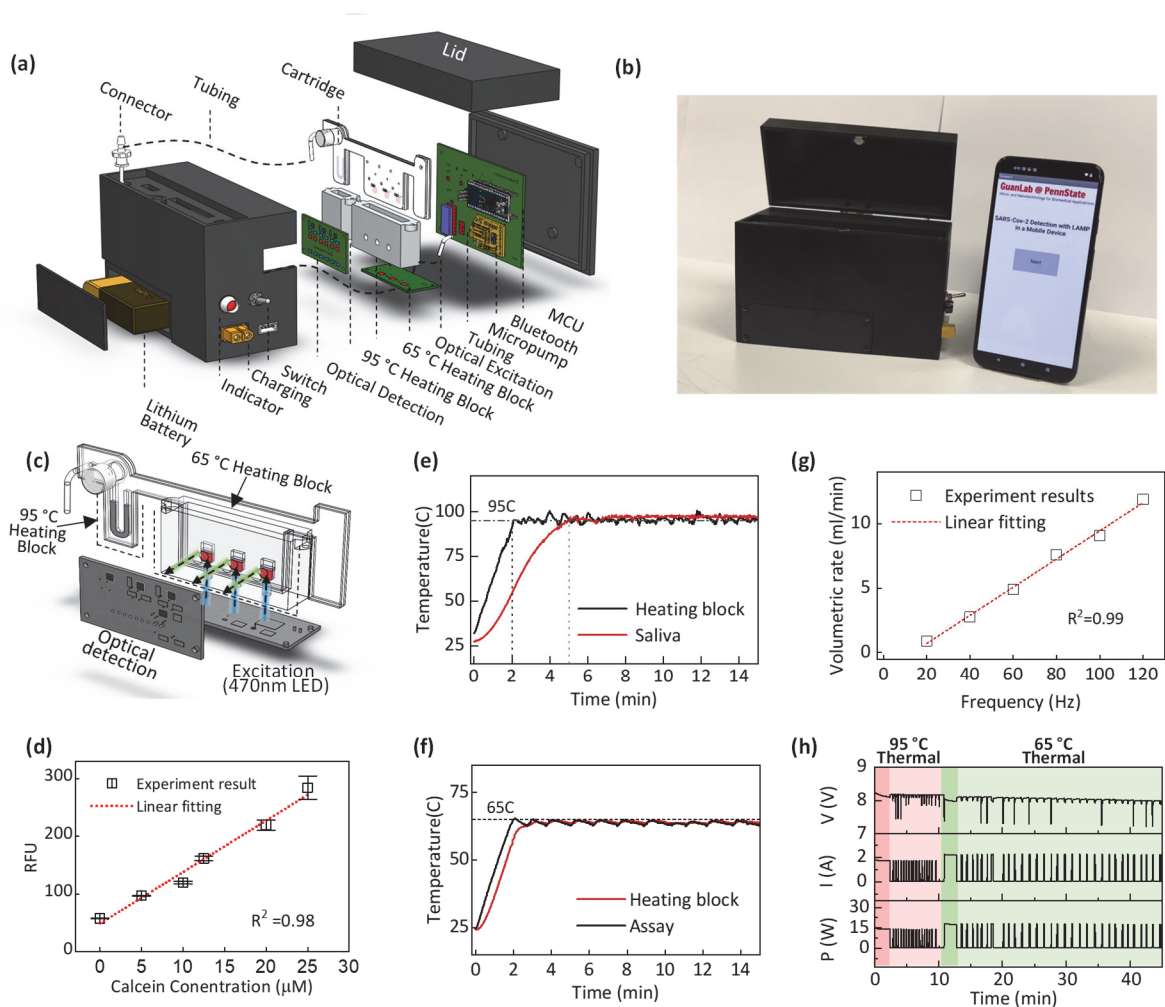


Figure 1. SLIDE Instrument design and validation. (a) Schematic of the SLIDE device showing components in an exploded view. The platform consists of five main functional modules: optical module (LED/optical sensor), thermal module (power resistor/thermal sensor), micro pump module, power supply module (battery), and data connectivity module (Bluetooth). Each module was controlled by a microcontroller on a customized PCB board. (b) Photograph of the SLIDE analyzer and the smartphone interface. (c) Schematic of the cartridge coordinated with the optical module and thermal module. (d) Characterization of the optical sensor using 40 μl calcein solution for 10 minutes of RFU recording. The optical sensor showed a linear response to the concentration of calcein from 0 to 25 μM . The temperature profile of the heating block and the liquid (saliva/assay) for (e) 95°C virus thermal lysis and (f) 65°C RT-LAMP reactions. (g) Characterization of the piezo pump frequency with the volumetric rate. (h) Power consumption characterization for one complete test.

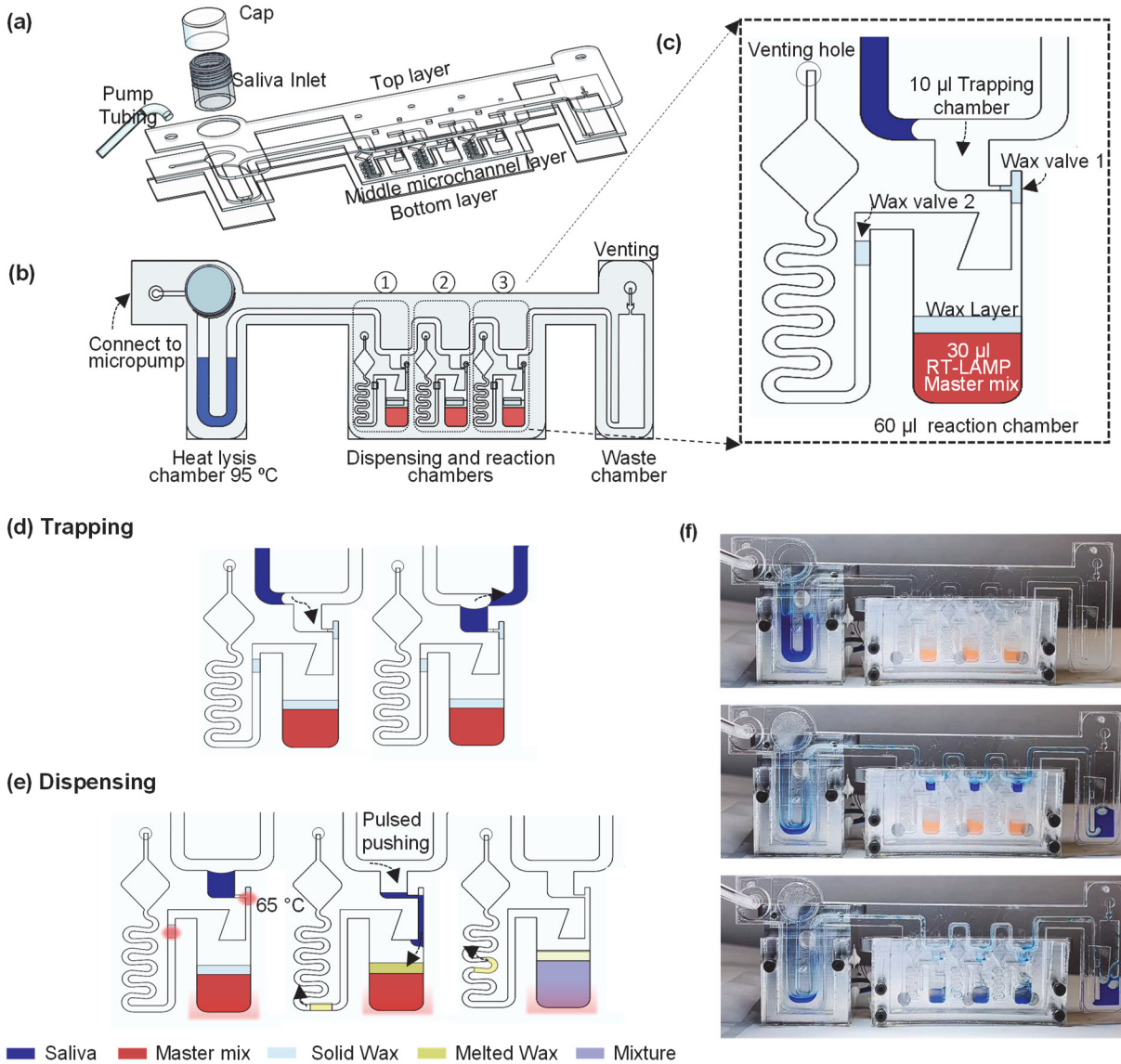


Figure 2. (a) Exploded view of the cartridge with three PMMA layers: top loading layer, middle microchannel layer, and bottom covering layer. (b) Assembled view of the cartridge includes a saliva collection chamber (250 µl), three dispensing and reaction chambers, and a waste chamber (300 µl). (c) One unit of the dispensing and reaction chamber comprise one trapping chamber (10 µl), two wax valves, one reaction chamber (60 µl) with preloaded RT-LAMP master mix and wax layer, as well as a venting hole for connect the atmosphere. Illustration of (d) trapping and (e) dispensing processes. (f) One example of sample trapping and dispensing processes (**Supplementary Video S2**). The blue liquid is the saliva mixed with the blue dye, and the orange liquid is the RT-LAMP master mixed with orange dye for better visualization.

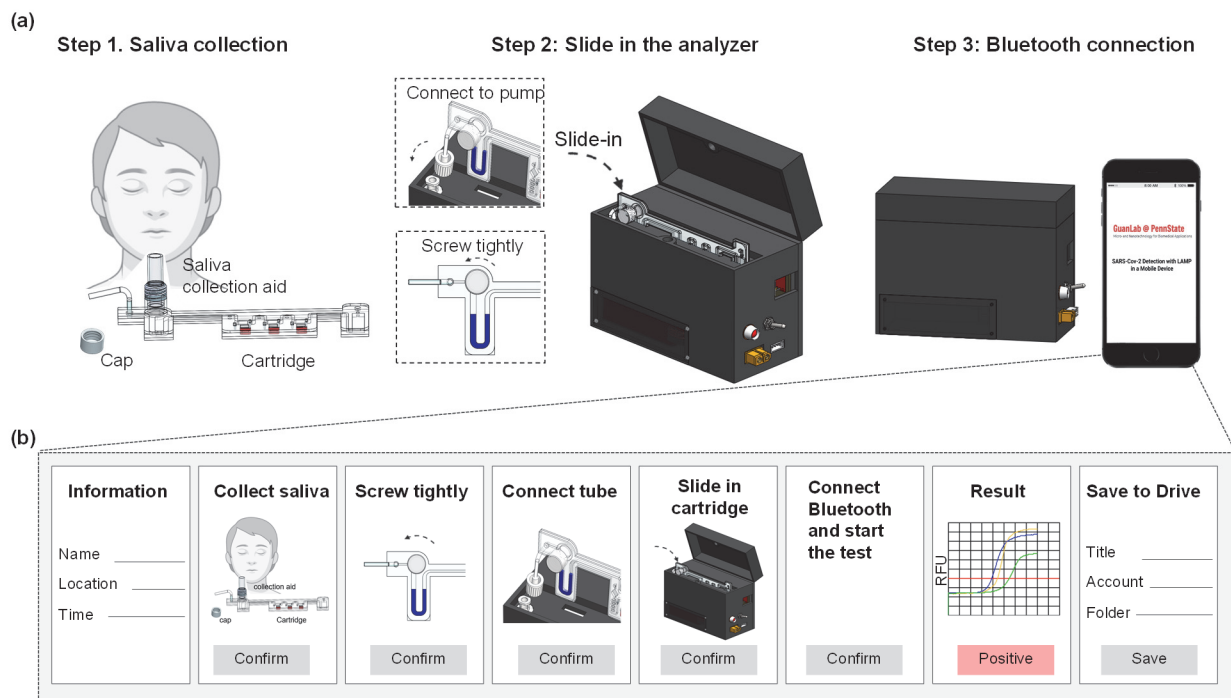


Figure 3. Overall SLIDE workflow. (a) Step 1: users self-collect $\sim 120 \mu\text{L}$ of saliva into a cartridge with the help of a saliva collection aid. Users tighten the screw cap and connect the Luer-lock to the micro pump. Step 2: Insert the cartridge into the analyzer and close the lid. Step 3: Connect the SLIDE analyzer with a smartphone through Bluetooth to initiate the test. (b) The step-by-step instruction of the APP interface includes personal information collection, sample collection guidance, Bluetooth connection, test initiation, data processing and communication.

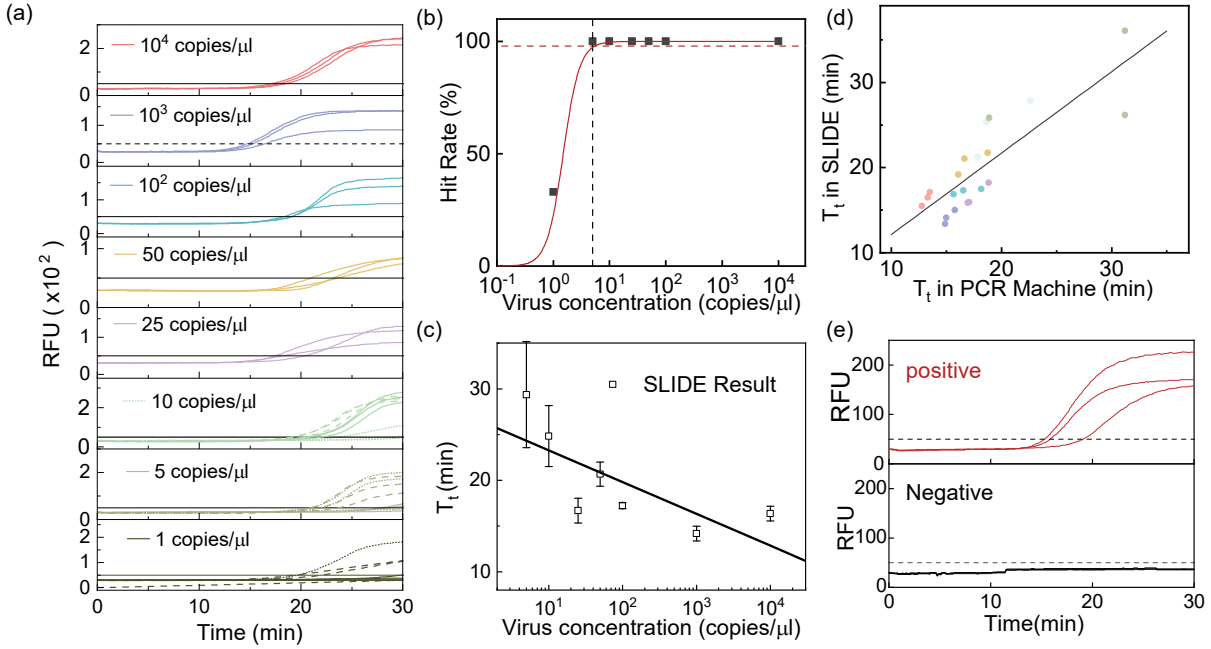


Figure 4. SLIDE device performance evaluation. (a) Real-time RT-LAMP results with different concentrations of the spiked saliva samples (1 virion/ μ l to 10^4 virion/ μ l) using a SLIDE analyzer. The threshold to classify an amplification curve as positive or negative was 50 RFU based on the NTC sample (**Supplementary Figure S5**). (b) The extracted hit rate at various virus particle concentrations to establish LoD. (c) The inversely proportional relationship between the threshold time (T_t) and virus particle concentration was obtained from the SLIDE analyzer. (d) The Pearson correlation analysis of the T_t between the manual operation with PCR machine and automatic method using SLIDE analyzer. (e) Two clinical samples, one known positive(top) and one known negative(bottom), were tested by the SLIDE device. The initial diagnosis was performed by the RT-PCR assay as the reference.

References

1. Harapan, H.; Itoh, N.; Yufika, A.; Winardi, W.; Keam, S.; Te, H.; Megawati, D.; Hayati, Z.; Wagner, A. L.; Mudatsir, M., Coronavirus disease 2019 (COVID-19): A literature review. *Journal of infection and public health* **2020**, *13* ((5)), 667-673.
2. Zhu, N.; Zhang, D.; Wang, W.; Li, X.; Yang, B.; Song, J.; Zhao, X.; Huang, B.; Shi, W.; Lu, R., A novel coronavirus from patients with pneumonia in China, 2019. *New England journal of medicine* **2020**, *382*, 727-733.
3. Kwok, H.; Briggs, K.; Tabard-Cossa, V., Nanopore Fabrication by Controlled Dielectric Breakdown. *PloS one* **2014**, *9* (3).
4. Chan, J. F.-W.; Yip, C. C.-Y.; To, K. K.-W.; Tang, T. H.-C.; Wong, S. C.-Y.; Leung, K.-H.; Fung, A. Y.-F.; Ng, A. C.-K.; Zou, Z.; Tsoi, H.-W., Improved molecular diagnosis of COVID-19 by the novel, highly sensitive and specific COVID-19-RdRp/Hel real-time reverse transcription-PCR assay validated in vitro and with clinical specimens. *Journal of clinical microbiology* **2020**, *58* (5), e00310-20.
5. Nagura-Ikeda, M.; Imai, K.; Tabata, S.; Miyoshi, K.; Murahara, N.; Mizuno, T.; Horiuchi, M.; Kato, K.; Imoto, Y.; Iwata, M., Clinical evaluation of self-collected saliva by quantitative reverse transcription-PCR (RT-qPCR), direct RT-qPCR, reverse transcription–loop-mediated isothermal amplification, and a rapid antigen test to diagnose COVID-19. *Journal of Clinical Microbiology* **2020**, *58* (9), e01438-20.
6. Tahamtan, A.; Ardebili, A., Real-time RT-PCR in COVID-19 detection: issues affecting the results. *Expert review of molecular diagnostics* **2020**, *20* (5), 453-454.
7. Vogels, C. B.; Brito, A. F.; Wyllie, A. L.; Fauver, J. R.; Ott, I. M.; Kalinich, C. C.; Petrone, M. E.; Casanovas-Massana, A.; Muenker, M. C.; Moore, A. J., Analytical sensitivity and efficiency comparisons of SARS-CoV-2 RT–qPCR primer–probe sets. *Nature microbiology* **2020**, *5* (10), 1299-1305.
8. Liao, W. T.; Hsu, M. Y.; Shen, C. F.; Hung, K. F.; Cheng, C. M., Home Sample Self - Collection for COVID - 19 Patients. *Advanced Biosystems* **2020**, *4* (11), 2000150.
9. Mak, G. C.; Cheng, P. K.; Lau, S. S.; Wong, K. K.; Lau, C.; Lam, E. T.; Chan, R. C.; Tsang, D. N., Evaluation of rapid antigen test for detection of SARS-CoV-2 virus. *Journal of Clinical Virology* **2020**, *129*, 104500.
10. Dinnes, J.; Deeks, J. J.; Berhane, S.; Taylor, M.; Adriano, A.; Davenport, C.; Dittrich, S.; Emperador, D.; Takwoingi, Y.; Cunningham, J., Rapid, point-of-care antigen and molecular-based tests for diagnosis of SARS-CoV-2 infection. *Cochrane Database of Systematic Reviews* **2021**, *Revision* (3).
11. Albert, E.; Torres, I.; Bueno, F.; Huntley, D.; Molla, E.; Fernández-Fuentes, M. Á.; Martínez, M.; Poujois, S.; Forqué, L.; Valdivia, A., Field evaluation of a rapid antigen test (Panbio™ COVID-19 Ag Rapid Test Device) for COVID-19 diagnosis in primary healthcare centres. *Clinical Microbiology and Infection* **2021**, *27* (3), 472. e7-472. e10.
12. Torres, I.; Poujois, S.; Albert, E.; Colomina, J.; Navarro, D., Evaluation of a rapid antigen test (Panbio™ COVID-19 Ag rapid test device) for SARS-CoV-2 detection in asymptomatic close contacts of COVID-19 patients. *Clinical Microbiology and Infection* **2021**, *27* (4), 636. e1-636. e4.
13. Scohy, A.; Anantharajah, A.; Bodéus, M.; Kabamba-Mukadi, B.; Verroken, A.; Rodriguez-Villalobos, H., Low performance of rapid antigen detection test as frontline testing for

COVID-19 diagnosis. *Journal of Clinical Virology* **2020**, *129*, 104455.

14. Brümmer, L. E.; Katzenschlager, S.; Gaeddert, M.; Erdmann, C.; Schmitz, S.; Bota, M.; Grilli, M.; Larmann, J.; Weigand, M. A.; Pollock, N. R., Accuracy of novel antigen rapid diagnostics for SARS-CoV-2: A living systematic review and meta-analysis. *PLoS medicine* **2021**, *18* (8), e1003735.
15. Taleghani, N.; Taghipour, F., Diagnosis of COVID-19 for controlling the pandemic: A review of the state-of-the-art. *Biosensors and Bioelectronics* **2021**, *174*, 112830.
16. Van Nguyen, H.; Nguyen, V. D.; Nguyen, H. Q.; Chau, T. H. T.; Lee, E. Y.; Seo, T. S., Nucleic acid diagnostics on the total integrated lab-on-a-disc for point-of-care testing. *Biosensors and Bioelectronics* **2019**, *141*, 111466.
17. Kojima, N.; Turner, F.; Slepnev, V.; Bacelar, A.; Deming, L.; Kodeboyina, S.; Klausner, J., Self-Collected Oral Fluid and Nasal Swabs Demonstrate Comparable Sensitivity to Clinician Collected Nasopharyngeal Swabs for Coronavirus Disease 2019 Detection. *Clinical Infectious Diseases* **2021**, *73* (9), e3106-e3109.
18. Williams, E.; Bond, K.; Zhang, B.; Putland, M.; Williamson, D. A., Saliva as a Noninvasive Specimen for Detection of SARS-CoV-2. *Journal of Clinical Microbiology* **2020**, *58* (8), e00776-20.
19. Pasomsub, E.; Watcharananan, S. P.; Boonyawat, K.; Janchompoo, P.; Wongtabtim, G.; Sukswan, W.; Sungkanuparph, S.; Phuphuakrat, A., Saliva sample as a non-invasive specimen for the diagnosis of coronavirus disease 2019: a cross-sectional study. *Clinical Microbiology and Infection* **2021**, *27* (2), 285. e1-285. e4.
20. Wyllie, A. L.; Fournier, J.; Casanovas-Massana, A.; Campbell, M.; Tokuyama, M.; Vijayakumar, P.; Warren, J. L.; Geng, B.; Muenker, M. C.; Moore, A. J., Saliva or nasopharyngeal swab specimens for detection of SARS-CoV-2. *New England Journal of Medicine* **2020**, *383* (13), 1283-1286.
21. Vogels, C. B.; Watkins, A. E.; Harden, C. A.; Brackney, D. E.; Shafer, J.; Wang, J.; Caraballo, C.; Kalinich, C. C.; Ott, I. M.; Fauver, J. R., SalivaDirect: A simplified and flexible platform to enhance SARS-CoV-2 testing capacity. *Med* **2021**, *2* (3), 263-280.
22. Ning, B.; Yu, T.; Zhang, S.; Huang, Z.; Tian, D.; Lin, Z.; Niu, A.; Golden, N.; Hensley, K.; Threton, B., A smartphone-read ultrasensitive and quantitative saliva test for COVID-19. *Science advances* **2021**, *7* (2), eabe3703.
23. Lalli, M. A.; Langmade, J. S.; Chen, X.; Fronick, C. C.; Sawyer, C. S.; Burcea, L. C.; Wilkinson, M. N.; Fulton, R. S.; Heinz, M.; Buchser, W. J., Rapid and extraction-free detection of SARS-CoV-2 from saliva by colorimetric reverse-transcription loop-mediated isothermal amplification. *Clinical chemistry* **2021**, *67* (2), 415-424.
24. Azzi, L.; Carcano, G.; Gianfagna, F.; Grossi, P.; Dalla Gasperina, D.; Genoni, A.; Fasano, M.; Sessa, F.; Tettamanti, L.; Carinci, F., Saliva is a reliable tool to detect SARS-CoV-2. *Journal of Infection* **2020**, *81* (1), e45-e50.
25. Butler-Laporte, G.; Lawandi, A.; Schiller, I.; Yao, M.; Dendukuri, N.; McDonald, E. G.; Lee, T. C., Comparison of saliva and nasopharyngeal swab nucleic acid amplification testing for detection of SARS-CoV-2: a systematic review and meta-analysis. *JAMA internal medicine* **2021**, *181* (3), 353-360.
26. Notomi, T.; Okayama, H.; Masubuchi, H.; Yonekawa, T.; Watanabe, K.; Amino, N.; Hase, T., Loop-mediated isothermal amplification of DNA. *Nucleic acids research* **2000**, *28* (12), e63-e63.
27. Parida, M.; Sannarangaiah, S.; Dash, P. K.; Rao, P.; Morita, K., Loop mediated isothermal

- amplification (LAMP): a new generation of innovative gene amplification technique; perspectives in clinical diagnosis of infectious diseases. *Reviews in medical virology* **2008**, *18* (6), 407-421.
28. Subramanian, S.; Gomez, R. D., An empirical approach for quantifying loop-mediated isothermal amplification (LAMP) using *Escherichia coli* as a model system. *PLoS One* **2014**, *9* (6), e100596.
 29. Roy, S.; Wei, S. X.; Ying, J. L. Z.; Safavieh, M.; Ahmed, M. U., A novel, sensitive and label-free loop-mediated isothermal amplification detection method for nucleic acids using luminophore dyes. *Biosensors and Bioelectronics* **2016**, *86*, 346-352.
 30. Tang, Z.; Choi, G.; Nouri, R.; Guan, W., Loop-Mediated Isothermal Amplification-Coupled Glass Nanopore Counting Toward Sensitive and Specific Nucleic Acid Testing. *Nano letters* **2019**, *19* (11), 7927-7934.
 31. Baek, Y. H.; Um, J.; Antigua, K. J. C.; Park, J.-H.; Kim, Y.; Oh, S.; Kim, Y.-i.; Choi, W.-S.; Kim, S. G.; Jeong, J. H., Development of a reverse transcription-loop-mediated isothermal amplification as a rapid early-detection method for novel SARS-CoV-2. *Emerging microbes & infections* **2020**, *9* (1), 998-1007.
 32. Ganguli, A.; Mostafa, A.; Berger, J.; Aydin, M. Y.; Sun, F.; de Ramirez, S. A. S.; Valera, E.; Cunningham, B. T.; King, W. P.; Bashir, R., Rapid isothermal amplification and portable detection system for SARS-CoV-2. *Proceedings of the National Academy of Sciences* **2020**, *117* (37), 22727-22735.
 33. Jiang, M.; Pan, W.; Arasthfer, A.; Fang, W.; Ling, L.; Fang, H.; Daneshnia, F.; Yu, J.; Liao, W.; Pei, H., Development and validation of a rapid, single-step reverse transcriptase loop-mediated isothermal amplification (RT-LAMP) system potentially to be used for reliable and high-throughput screening of COVID-19. *Frontiers in cellular and infection microbiology* **2020**, *10*, 331.
 34. Kitagawa, Y.; Orihara, Y.; Kawamura, R.; Imai, K.; Sakai, J.; Tarumoto, N.; Matsuoka, M.; Takeuchi, S.; Maesaki, S.; Maeda, T., Evaluation of rapid diagnosis of novel coronavirus disease (COVID-19) using loop-mediated isothermal amplification. *Journal of Clinical Virology* **2020**, *129*, 104446.
 35. Lamb, L. E.; Bartolone, S. N.; Ward, E.; Chancellor, M. B., Rapid detection of novel coronavirus (COVID19) by reverse transcription-loop-mediated isothermal amplification. Available at SSRN 3539654 **2020**, <http://dx.doi.org/10.2139/ssrn.3539654>.
 36. Lu, R.; Wu, X.; Wan, Z.; Li, Y.; Zuo, L.; Qin, J.; Jin, X.; Zhang, C., Development of a novel reverse transcription loop-mediated isothermal amplification method for rapid detection of SARS-CoV-2. *Virologica Sinica* **2020**, *35* (3), 344-347.
 37. Park, G.-S.; Ku, K.; Baek, S.-H.; Kim, S.-J.; Kim, S. I.; Kim, B.-T.; Maeng, J.-S., Development of reverse transcription loop-mediated isothermal amplification assays targeting severe acute respiratory syndrome coronavirus 2 (SARS-CoV-2). *The Journal of Molecular Diagnostics* **2020**, *22* (6), 729-735.
 38. Yan, C.; Cui, J.; Huang, L.; Du, B.; Chen, L.; Xue, G.; Li, S.; Zhang, W.; Zhao, L.; Sun, Y., Rapid and visual detection of 2019 novel coronavirus (SARS-CoV-2) by a reverse transcription loop-mediated isothermal amplification assay. *Clinical Microbiology and Infection* **2020**, *26* (6), 773-779.
 39. Yu, L.; Wu, S.; Hao, X.; Dong, X.; Mao, L.; Pelechano, V.; Chen, W.-H.; Yin, X., Rapid detection of COVID-19 coronavirus using a reverse transcriptional loop-mediated isothermal amplification (RT-LAMP) diagnostic platform. *Clinical chemistry* **2020**, *66* (7), 975-977.

40. Zhu, X.; Wang, X.; Han, L.; Chen, T.; Wang, L.; Li, H.; Li, S.; He, L.; Fu, X.; Chen, S., Multiplex reverse transcription loop-mediated isothermal amplification combined with nanoparticle-based lateral flow biosensor for the diagnosis of COVID-19. *Biosensors and Bioelectronics* **2020**, *166*, 112437.
41. Tang, Z.; Nouri, R.; Dong, M.; Yang, J.; Greene, W.; Zhu, Y.; Yon, M.; Nair, M. S.; Kuchipudi, S. V.; Guan, W., Rapid detection of novel coronavirus SARS-CoV-2 by RT-LAMP coupled solid-state nanopores. *Biosensors and Bioelectronics* **2022**, *197*, 113759.
42. Liu, T.; Choi, G.; Tang, Z.; Kshirsagar, A.; Politza, A. J.; Guan, W., Fingerpick Blood-Based Nucleic Acid Testing on A USB Interfaced Device towards HIV self-testing. *Biosensors & bioelectronics* **2022**, *209*, 114255.
43. Choi, G.; Prince, T.; Miao, J.; Cui, L.; Guan, W., Sample-to-answer palm-sized nucleic acid testing device towards low-cost malaria mass screening. *Biosensors and Bioelectronics* **2018**, *115*, 83-90.
44. Holstein, C. A.; Griffin, M.; Hong, J.; Sampson, P. D., Statistical method for determining and comparing limits of detection of bioassays. *Analytical chemistry* **2015**, *87* (19), 9795-9801.

For TOC only

